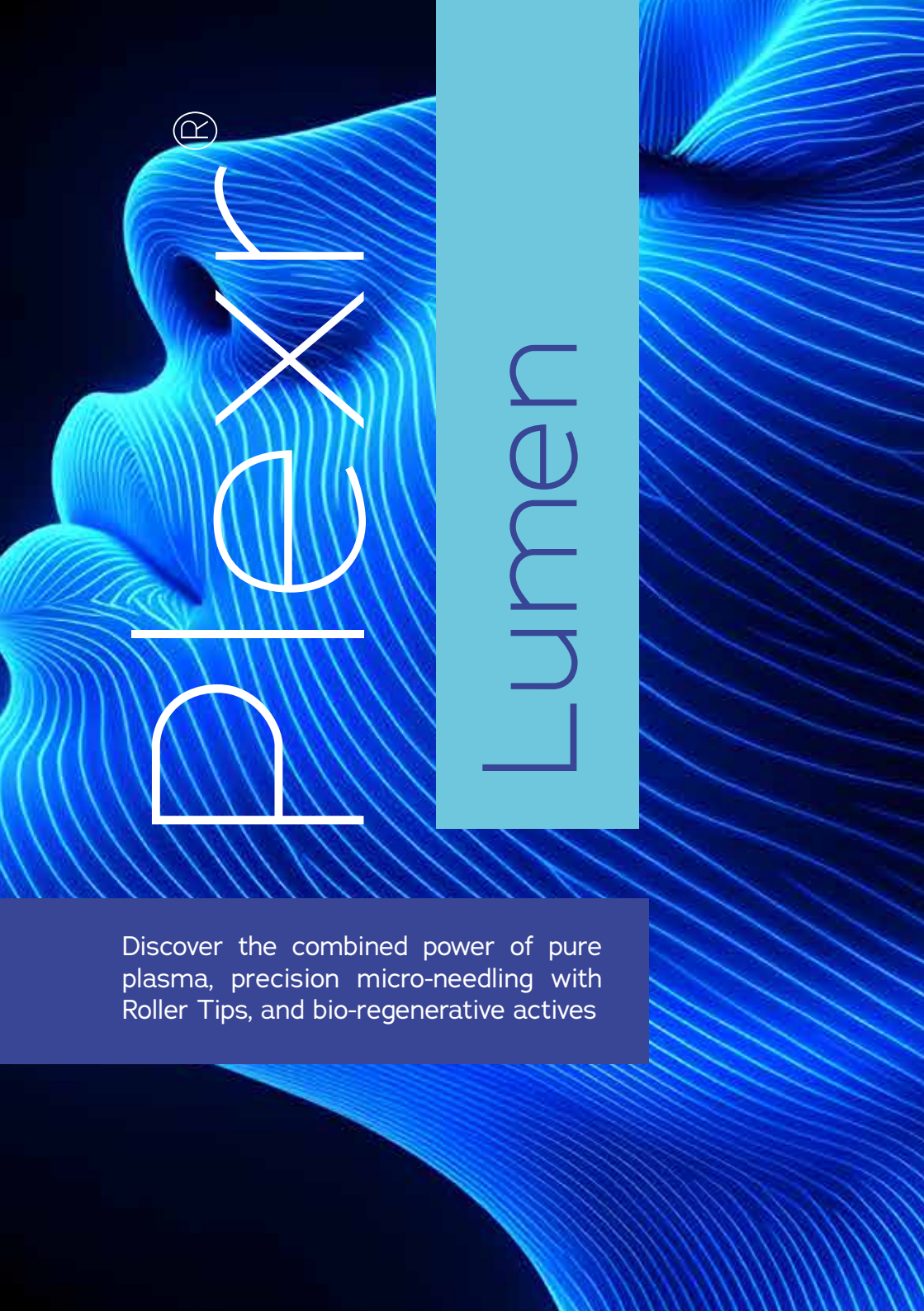

plexr®
Lumen
Discover the combined power of pure plasma, precision micro-needling with Roller Tips, and bio-regenerative actives

WHAT IS PLEXR® LUMEN

An innovative treatment that combines Plexr® energy with precision micro-needling to activate regenerative skin processes and enhance the delivery of functional bio-active ingredients



EXPECTED RESULTS

- Brighter, more even skin
- Visible reduction of fine lines and enlarged pores
- Increased elasticity and hydration

Experience the effectiveness of a protocol that merges micro-plasma, micro-needling, and bio-regenerative actives for long-lasting, scientifically-supported rejuvenation



HOW IT WORKS

During micro-needling with Roller Tips, Plexr® energy is delivered directly through the micro-needles, combining controlled energetic stimulation with mechanical micro-channeling.
This dual action increases permeability, enhances regenerative signaling, and promotes collagen reorganization.

Bio-regenerative actives — including Plexsomes — reach deeper epidermal layers, integrating with the known biological effects of Plexr®.

The result is a comprehensive regenerative process: preparation, stimulation, and infusion working in harmony for an advanced minimally invasive treatment.

THE SYSTEM

ROLLER HOLDER

The engineered connector linking Roller Tips to Plexr®, ensuring secure coupling and optimal energy transmission — transforming the roller into an active therapeutic tool

Compatible with ALL Plexr® models

ROLLER TIPS

Gold-coated stainless-steel micro-needles designed for precise, uniform micro-channeling and efficient delivery of bio-active ingredients.

Lightweight, sterile, single-use, and available in multiple depths (0.25 / 0.5 / 1.0 mm).

Roller Tips can be used in two ways:

- with Plexr® energy alone, for mechanical micro-channeling combined with controlled energetic stimulation
- or together with Plexr® energy and Plexsomes (or other active bio-regeneratives), when a deeper and more complete regenerative biological response is desired

TREATMENT PROTOCOL

STEP 1

- Cleanse the skin and disinfect with a non-alcoholic antiseptic.
- (Optional) Perform a light plasma peeling with BioTips to improve turnover and uniformity.

STEP 2

- Insert the Roller Holder and apply the selected Roller Tip.
- Perform an initial light and even pass.
- Apply Plexsomes or another regenerative active.
- Repeat vertical, horizontal, and diagonal passes to maximize transdermal diffusion.
- Add additional product as needed to maintain hydration during the treatment (recommended: approx. 2–2.5 mL for full face).

STEP 3

- Apply a soothing cream (e.g., zinc oxide), optionally delivered with BioTips, to stabilize the skin surface.

Treatment Plan

- 3 in-clinic sessions (Day 0, Day 15, Day 30) with an optional 4th session (Day 45)
- 2 at-home supportive treatments (scrub and mask) between sessions
- Maintenance recommended with a touch-up session after 6 months

Note

- Select the Roller Tip depth and Plexr® handpiece based on: skin type, sensitivity, and treatment goal.
- Topical anesthetic may be used 30 minutes before treatment if the patient is sensitive.

BENEFITS & WHY CHOOSE IT

Stimulates collagen production and improves overall skin quality

Enhances the transdermal delivery of exosomes, peptides, and functional serums

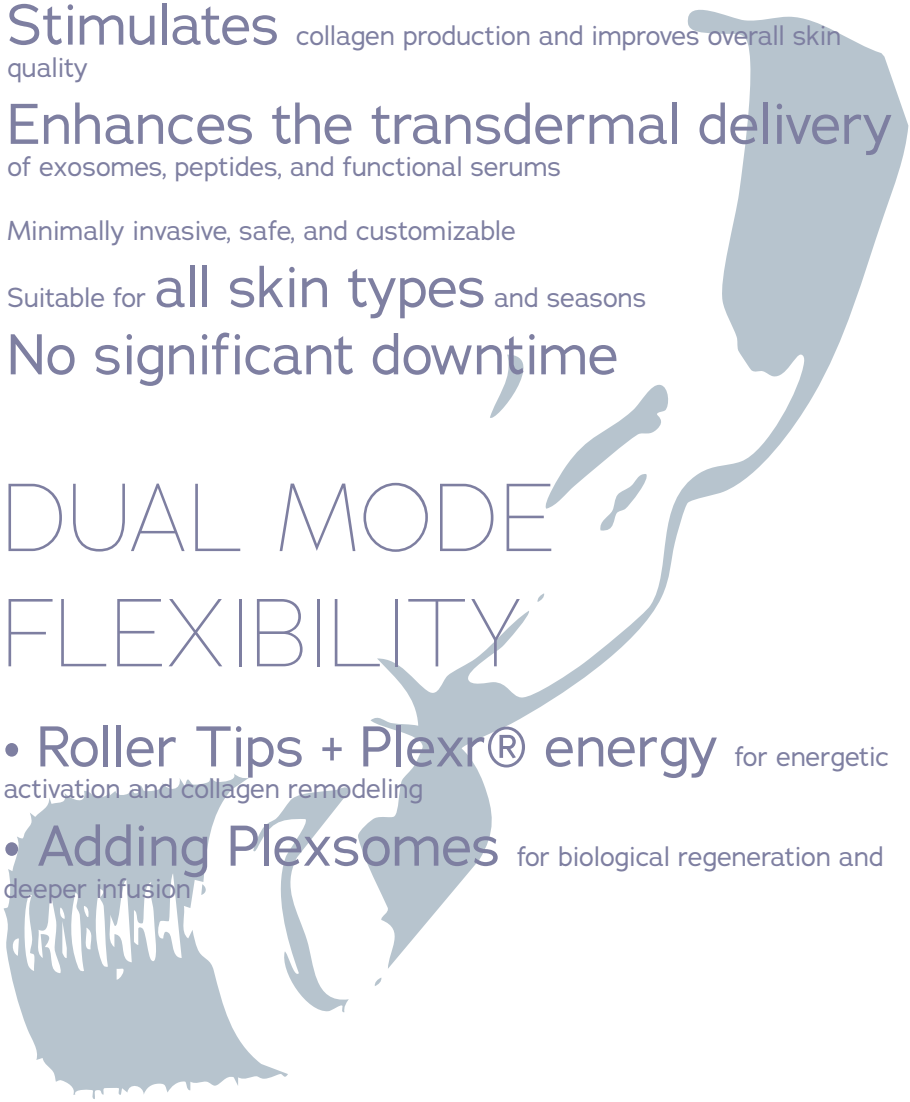
Minimally invasive, safe, and customizable

Suitable for **all skin types** and seasons

No significant downtime

DUAL MODE FLEXIBILITY

- **Roller Tips + Plexr® energy** for energetic activation and collagen remodeling
- **Adding Plexsomes** for biological regeneration and deeper infusion





www.gruppogmv.com